




# 1 How to log in the system



1. Press the Log in button

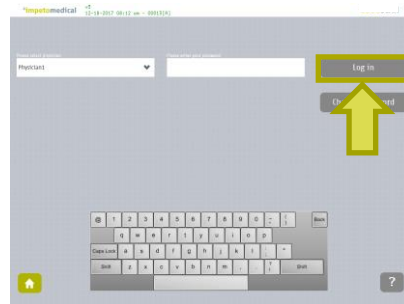

2. Select a user and enter the password. Press Log in.




3. The user is connected


# 2 How to start a scan for a new patient



4. Press the New Patient button.

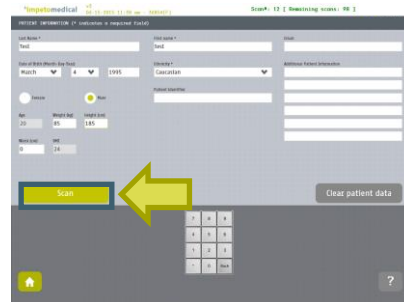

5. Fill in demographics data and press Scan.

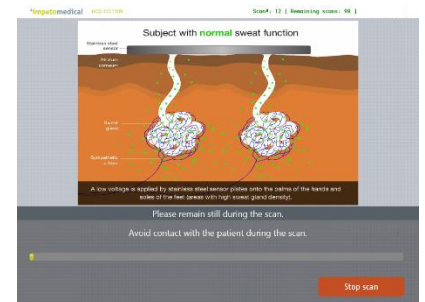

6. Scan in progress...


# 3 How to view the results



7. Results: Press Physician report or Patient report to display a report.



9. Patient report ("Optional")

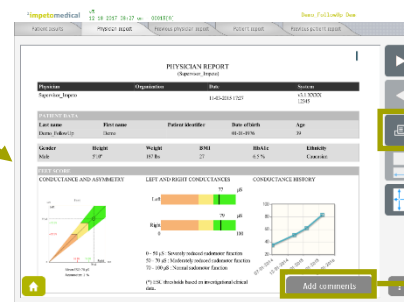

8. Physician report

Print

Add comments

Print


Add comments




## 4 How to retrieve data for a patient




10. Press Patient History

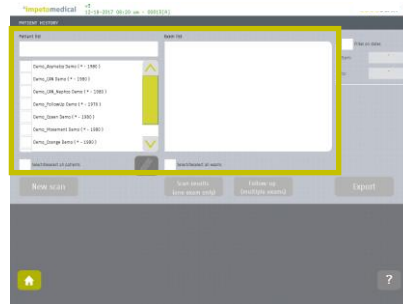


11. Select the patient and select one or several exams

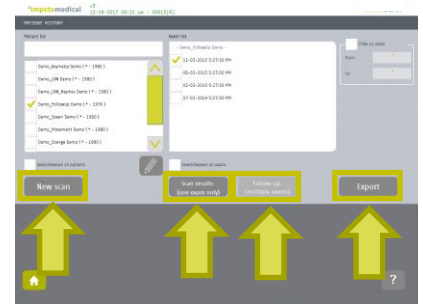

12. Select an action

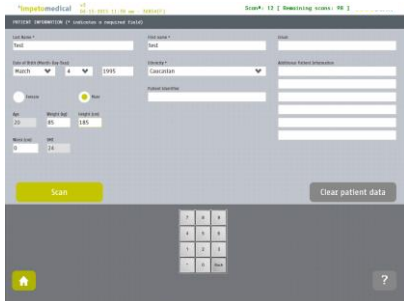


13. NEW SCAN: Start a new scan for the selected patient (only one patient selected)

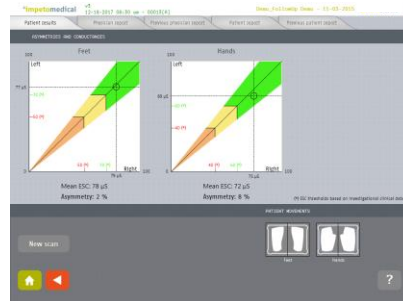

14. SCAN RESULTS: Display the results (only one exam selected)



15. FOLLOW-UP: Display the follow-up chart for the selected patient (Two or more exams selected for one patient)

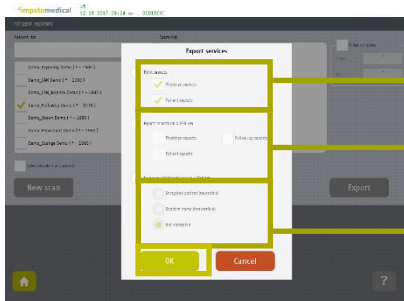


16. EXPORT : Display the export functions


Print reports (physician and/or patient)

Export reports (physician, patient, and follow-up) on a USB key (cf. \ExportReports directory)


Compile a CSV file for selected exams and export it to a USB key (cf. \CompiledCSV directory).  
These files can be opened with Excel or Open Office.

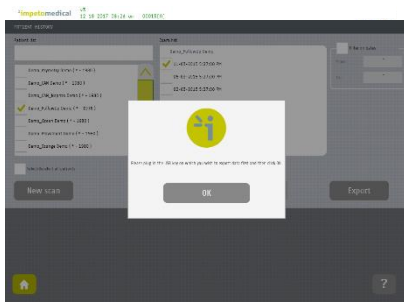

17. Plug in a USB key and validate

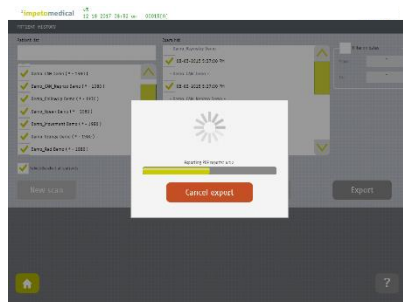

18. The files are being exported on the USB key

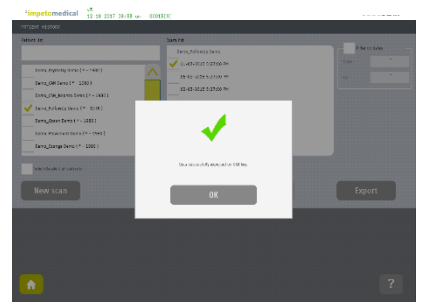

19. The files are successfully exported